

CureMed Nordic AB is a Swedish medtech company with head office and warehouse in Gothenburg. In our range, you will find world-leading products within healthcare and medical technology, many of which can be customized in production. Our goal is to always offer innovative products from leading manufacturers in a cost-effective manner to the healthcare sector.

Replacement Clamp™

- Replacement clamp is a repair clamp for dialysis catheters*.
- Replacement clamp makes it possible to replace a broken clamp on the dialysis catheter without removing the luer lock connector.
- New ergonomic design.
- Manufactured in Sweden.
- CE certified.




Step 1:

Remove the old clamp. Slide in the Replacement clamp from the side.

Step 2:

Fix the tube in the notch of the Replacement clamp. The clamp is now ready to use.

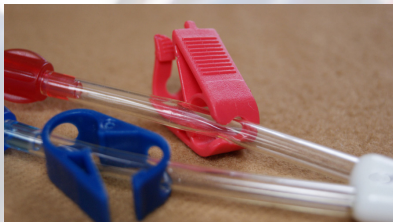

**Replacement clamp is for temporary use only. A damaged catheter should always be replaced. Follow the instructions of your catheter supplier carefully.*

ORDERING INFORMATION

Product code	Description	Packaging
RCR01	Replacement clamp, red	2 pcs / pouch
RCB02	Replacement clamp, blue	2 pcs / pouch